

Airport Terminal Wheelchair

CHECK IN TO AIRCRAFT DOOR

The Melrose Airport Terminal Wheelchair transports wheelchair passengers from the check-in to the aircraft door. An assistant is required to push this wheelchair so it has a full width padded push bar for comfortable one or two hand pushing.

Fixed frame (non folding) wheelchairs provide longer life cycles because the frames have no moving parts. Frames are hand made from aircraft grade aluminium and is lightweight and rust resistant. Many of our components are made inhouse to ensure we maintain exact tolerances from the highest quality materials. All four wheels are solid, these provide less rolling resistance for easier pushing and maneuverability, they also require very little maintenance over their life.

Cushions are required for skin protection and comfort, the skin on a disabled passenger can be susceptible to pressure areas, abrasion and bruising, even after a short period of time on a hard surface. Armrests provide support for passengers with poor balance, they also allow the passenger to pressure relief. Armrests lift up to allow transfers from the side. A Hip Seat-belt provides support and extra security for passengers with poor balance.

STANDARD SPECIFICATIONS


- Lightweight Aluminum Frame
- Seat Width 20.5", Depth 17"
- 1" Cushion with wipe to clean Cover
- Adjustable height flip up Foot-Plates
- Lift-back Armrests
- 6" solid wide Castors
- 18" solid rear Wheels
- Hand operated Drum Brakes
- Black fabric Sling and Back
- Seat-belt with Car style Buckle
- Under-chair Luggage Carrier
- Wide range of Colours
- Max user weight 120kg
- Five-year Frame Warranty

OPTIONS AND ACCESSORIES

- Titanium frame
- Fluro and Metallic colours
- Seat-belt with Velcro Fastening
- Seat-belt with Plastic Clasp Fastening
- Airline or Airport logo printed on Seat Back



Lift-back Armrests

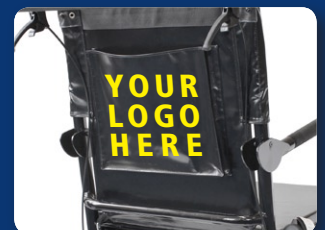

Logo on Seat Back